

Compact
Powerful
Precise



ColorFlex[®] L2 Tomato

*Level UP your Color Confidence with the new
ColorFlex L2 Tomato Spectrophotometer.*

As simple as you want, as powerful as you need.

Elevate your color accuracy and user experience like never before. Compact and ergonomic build requires minimal space while offering enhanced durability and reliability. Designed to streamline and expedite your quality control processes, ColorFlex L2 Tomato Spectrophotometer is expertly crafted to measure the color of tomatoes in all their processed forms, including paste/purée, sauce, ketchup, juice, and fresh tomatoes.




Explore the Next-Gen ColorFlex®...L2 Tomato



#1 Color Measurement Solution for Tomatoes

HunterLab Tomato Reference Standard: HunterLab, the creators of the original tomato color reference standard, offers the only reference tile that traces back to the discontinued BCR and USDA standards, ensuring the most accurate red color measurements.

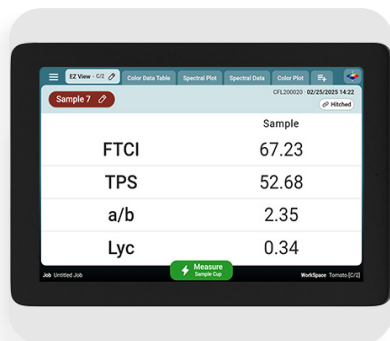


All the Tomato Scores and Indices You Need

Measure Tomato in Every Form: Tomato Paste Score (TPS), Tomato Sauce Score (TSS), Tomato Catsup Score (TCS), Tomato Juice Score (TJS), a/b ratio, Fresh Tomato Color Index (FTCI) & Lycopene index.

45°a:0° Annular Illuminated Optical Geometry

Superior Consistency and Color Accuracy: Enhancing the efficiency of Quality Control workflows and delivery of superior quality products repeatedly.

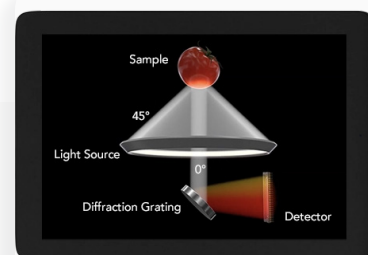

Compact with Sealed Spillproof Case

Space-Saving and Easy Maintenance: Occupies minimal space and is designed to withstand the demands of busy lab environments.



Tomato Paste
Tomatoes

Tomato Sauce
TOMATO KETCHUP

tomato juice
Tomato Soup



Powerful Standalone Functionality

No computer required: With huge onboard processing and storage capacity, ColorFlex L2 Tomato can function as an independent color workstation. Connect monitor, keyboard, and mouse if desired.




Modern Tablet User Interface with Touchscreen Color Display

Durable and Efficient: Facilitates accurate color measurement in seconds, with a robust industrial glass design suitable for harsh environments.




New & Improved Color Software

Essentials 2.0: Take precise and reliable color measurement in seconds with next-gen user-friendliness and seamlessly analyze multiple tristimulus and spectral data sets. It's that powerful.

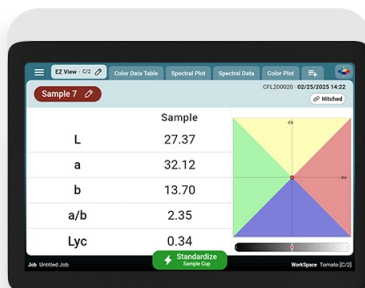

Self-Guided and Intuitive Out-of-Box Startup

Time-Saving and User-Friendly: Ready to use straight out of the box with an integrated wizard for quick training, reducing operator time and easing concerns related to staff training and turnover.

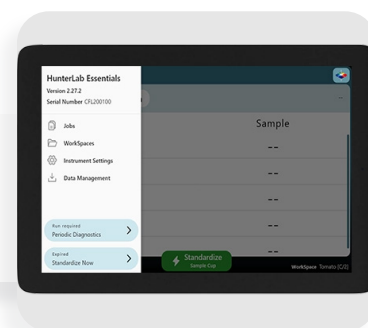

Built-in Camera for Sample Viewing and Screen Recording

Precision and Communication: Ensures accurate sample positioning to eliminate measurement discrepancies and allows easy sharing of recordings for swift quality control decisions.



Tomato Paste
Tomatoes


Tomato Sauce
TOMATO KETCHUP

tomato juice
Tomato Soup



Built-in Advanced Diagnostics

On-Board Diagnostic Health Check: Ensures the HunterLab ColorFlex L2 Tomato is always operating at peak performance, continuing to build on the legendary history of reliability of the ColorFlex family of instruments.

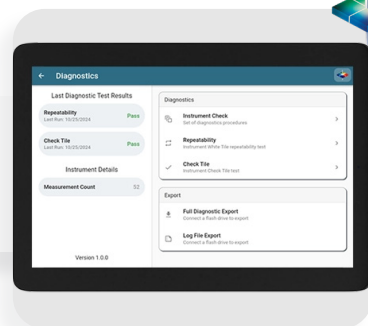

Advanced Data Export and Communication

Full suite of data communication options: Features USB ports, HDMI, and support for mouse and keyboard, enabling easy data export/printing.

Connectivity with LIMS and SPC systems

Simplifies QC workflows: Aids in trend analysis and workflow simplification across manufacturing sites.



CONTACT US

We've worked with thousands of companies, but you'll feel like our one and only. If color accuracy is crucial to your business, you've come to the right place—with the right people. Let's get in touch. Ask a question, request a quote or just say hello!

OUR OFFICES & DISTRIBUTION CENTERS

HunterLab

11491 Sunset Hills Road
Reston, VA 20190-5280
Website: www.hunterlab.com
Email: sales@hunterlab.com
Ph: 703-471-6870
Fax: 703-471-4237

Corporate Headquarters

HunterLab Europe

Christian Jansen, Managing Director
Dr.-August-Einsele-Ring 15
82418 Murnau, Germany
Website: www.hunterlab.com
Email: Christian.Jansen@hunterlab.com
Ph: +49 8841 9464

European Office

HunterLab Asia-Pacific

Heng Su Kiang, APAC Sales Manager
28-2, Plaza KLTS
No. 99, Jalan Gombak
Setapak
53000 Wilayah Persekutuan Kuala Lumpur
Malaysia
Website: www.hunterlab.com
Email: sales-asiapacific@hunterlab.com
Ph: +(60) 340 324 749

HunterLab Asia-Pacific

Leo Cheung, Regional Manager
Flat E, 13F Leahander Centre
28 Wang Wo Tsai Street Tsuen Wan
New Territories Hong Kong
Website: www.hunterlab.com
Email: sales-asiapacific@hunterlab.com
Ph: +852-2406 1723

Asia-Pacific Offices

For more information, please contact HunterLab at 703-471-6870, sales@hunterlab.com or visit www.hunterlab.com